# An R package to process LC/MS metabolomic data: MAIT (Metabolite Automatic Identification Toolkit)

Francesc Fernández-Albert, Rafael Llorach, Cristina Andrés-Lacueva, Alexandre Perera

May 1, 2024

# 1 Abstract

Processing metabolomic liquid chromatography and mass spectrometry (LC/MS) data files is time consuming. Currently available R tools allow for only a limited number of processing steps and online tools are hard to use in a programmable fashion. This paper introduces the metabolite automatic identification toolkit MAIT package, which allows users to perform end-to-end LC/MS metabolomic data analysis. The package is especially focused on improving the peak annotation stage and provides tools to validate the statistical results of the analysis. This validation stage consists of a repeated random sub-sampling cross-validation procedure evaluated through the classification ratio of the sample files. MAIT also includes functions that create a set of tables and plots, such as principal component analysis (PCA) score plots, cluster heat maps or boxplots. To identify which metabolites are related to statistically significant features, MAIT includes a metabolite database for a metabolite identification stage.

# 2 Introduction

Liquid Chromatography and Mass Spectrometry (LC/MS) is an analytical instrument widely used in metabolomics to detect molecules in biological samples. It breaks the molecules down into pieces, some of which are detected as peaks in the mass spectrometer. Metabolic profiling of LC/MS samples basically consists of a peak detection and signal normalisation step, followed by multivariate statistical analysis such as principal components analysis (PCA) and univariate statistical tests such as ANOVA .

As analysing metabolomic data is time consuming, a wide array of software tools are available, including commercial tools such as Analyst® software. There are programmatic R packages, such as XCMS to detect peaks or CAMERA package and AStream, which cover only peak annotation. Another category of free tools available consists of those having online access through a graphical user interface (GUI), such as XCMS Online (http://xcmsonline.scripps.edu) or MetaboAnalyst, both extensively used.

These online tools are difficult to use in a programmable fashion. They are also designed and programmed to be used step by step with user intervention, making it difficult to set

up metabolomic data analysis workflow. These R packages involve only a part of the entire metabolomic analysis process. Although there are specific R packages whose objective is peak annotation, this is still an issue in analysing LC/MS metabolomic data.

We introduce a new R package called metabolite automatic identification toolkit (MAIT) for automatic LC/MS analysis. The goal of the MAIT package is to provide an array of tools for programmable metabolomic end-to-end analysis. It consequently has special functions to improve peak annotation through the processes called biotransformations. Specifically, MAIT is designed to look for statistically significant metabolites that separate the classes in the data.

# 3 Methodology

The main processing steps for metabolomic LC/MS data include the following stages: peak detection, peak annotation and statistical analysis. In the peak detection stage, the objective is to detect the peaks in the LC/MS sample files. The peak annotation stage identifies the metabolites in the metabolomic samples better by increasing the chemical and biological information in the data set. A statistical analysis step is essential to obtain significant sample features. All these 3 steps are covered in the MAIT workflow.

## 3.1 Peak Detection

Peak detection in metabolomic LC/MS data sets is a complex issue for which several approaches have been developed. Two of the most well established techniques are matched filter and the centWave algorithm. MAIT can use both algorithms through the XCMS package.

## 3.2 Peak Annotation

The MAIT package uses 3 complementary steps in the peak annotation stage.

- The first annotation step uses a peak correlation distance approach and a retention time window to ascertain which peaks come from the same source metabolite, following the procedure defined in CAMERA package. The peaks within each peak group are annotated following a reference adduct/fragment table and a mass tolerance window.
- The second step uses a mass tolerance window inside the peak groups detected in the first step to look for more specific mass losses called biotransformations. To do this, MAIT uses a predefined biotransformation table where the biotransformations we want to find are saved. A user-defined biotransformation table can be set as an input following the procedure defined in Section (4.6).
- The third annotation step is the metabolite identification stage, in which a predefined metabolite database is mined to search for the significant masses, also using a tolerance window. This database is the Human Metabolome Database (HMDB), 2009/07 version.

## 3.3 Statistical Analysis

The objective of analysing metabolomic profiling data is to obtain the statistically significant features that contain the highest amount of class-related information. To gather these features, MAIT applies standard univariate statistical tests (ANOVA or Student's t-test) to every feature and selects the significant set of features by setting up a user-defined threshold P-value. Bonferroni multiple test correction can be applied to the resulting P-values. We propose a validation test to quantify how well the data classes are separated by the statistically significant features. The separation is validated through a repeated random sub-sampling cross-validation using partial least squares and discriminant analysis (PLS-DA), support vector machine (SVM) with a radial Kernel and K-nearest neighbours (KNN). Overall and class-related classification ratios are obtained to evaluate the class-related information of the significant features.

# 4 Using MAIT

The data files for this example are a subset of the data used in reference, which are freely distributed through the XCMS package. In these data there are 2 classes of mice: a group where the fatty acid amide hydrolase gene has been suppressed (class knockout or KO) and a group of wild type mice (class wild type or WT). There are 6 spinal cord samples in each class. In the following, the MAIT package will be used to read and analyse these samples using the main functions discussed in Section ??. The significant features related to each class will be found using statistical tests and analysed through the different plots that MAIT produces.

# 4.1 Data Import

Each sample class file should be placed in a directory with the class name. All the class folders should be placed under a directory containing only the folders with the files to be analysed. In this case, 2 classes are present in the data. An example of correct file distribution using the example data files is shown in Figure 1.

#### 4.2 Peak Detection

Once the data is placed in 2 subdirectories of a single folder, the function sampleProcessing() is run to detect the peaks, group the peaks across samples, perform the retention time correction and carry out the peak filling process. As function sampleProcessing() uses the XCMS package to perform these 4 processing steps, this function exposes XCMS parameters that might be modified to improve the peak detection step. A project name should be defined because all the tables and plots will be saved in a folder using that name. For example, typing project = "project\_Test", the output result folder will be "Results\_project\_Test".

By choosing "MAIT\_Demo" as the project name, the peak detection stage can be launched by typing:

> library(MAIT)

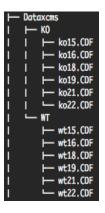


Figure 1: Example of the correct sample distribution for MAIT package use. Each sample file has to be saved under a folder with its class name.

```
> library(faahK0)
> cdfFiles<-system.file("cdf", package="faahK0", mustWork=TRUE)
> MAIT <- sampleProcessing(dataDir = cdfFiles, project = "MAIT_Demo",
+ snThres=2,rtStep=0.03)

Peak detection done
Retention time correction done
Peak grouping after samples done
Missing Peak integration done</pre>
```

After having launched the sampleProcessing function, peaks are detected, they are grouped across samples and their retention time values are corrected. A short summary in the R session can be retrieved by typing the name of the MAIT-class object.

#### > MAIT

```
A MAIT object built of 12 samples
The object contains 6 samples of class KO
The object contains 6 samples of class WT
```

The result is a MAIT-class object that contains information about the peaks detected, their class names and how many files each class contains. A longer summary of the data is retrieved by performing a summary of a MAIT-class object. In this longer summary version, further information related to the input parameters of the whole analysis is displayed. This functionality is especially useful in terms of traceability of the analysis.

```
> summary(MAIT)
```

```
A MAIT object built of 12 samples
The object contains 6 samples of class KO
The object contains 6 samples of class WT
```

```
Parameters of the analysis:
dataDir
                     "/home/biocbuild/bbs-3.20-bioc/R/site-library/faahKO/cdf"
snThres
Sigma
                     "2.12332257516562"
mzSlices
                     "0.3"
retcorrMethod
                     "loess"
                     "density"
groupMethod
bwGroup
                     "3"
                     "0.25"
mzWidGroup
filterMethod
                     "centWave"
rtStep
                     "0.03"
                     "0"
nSlaves
project
                     "MAIT_Demo"
ppm
                     "10"
                     "0.5"
minfrac
fwhm
                     "30"
family1
                     "gaussian"
                     "symmetric"
family2
                     "0.2"
span
centWave peakwidth1 "5"
centWave peakwidth2 "20"
```

#### 4.3 Peak Annotation

The next step in the data processing is the first peak annotation step, which is performed through the peakAnnotation(). If the input parameter adductTable is not set, then the default MAIT table for positive polarisation will be selected. However, if the adductTable parameter is set to "negAdducts", the default MAIT table for negative fragments will be chosen instead. peakAnnotation function also creates an output table (see Table ??) containing the peak mass (in charge/mass units), the retention time (in minutes) and the spectral ID number for all the peaks detected. A call of the function peakAnnotation may be:

```
> MAIT <- peakAnnotation(MAIT.object = MAIT,corrWithSamp = 0.7,
+ corrBetSamp = 0.75,perfwhm = 0.6)</pre>
```

WARNING: No input adduct/fragment table was given. Selecting default MAIT table for posit Set adductTable equal to negAdducts to use the default MAIT table for negative polarity Start grouping after retention time.

Created 332 pseudospectra.

Spectrum build after retention time done

Generating peak matrix!

```
Run isotope peak annotation
% finished: 10 20 30
                       40 50 60 70 80 90 100
Found isotopes: 104
Isotope annotation done
Start grouping after correlation.
Generating EIC's ..
Calculating peak correlations in 332 Groups...
 % finished: 10 20 30
                       40
                          50
                               60 70 80 90 100
Calculating peak correlations across samples.
                               60 70 80
 % finished: 10 20 30 40 50
                                          90
                                              100
Calculating isotope assignments in 332 Groups...
 % finished: 10 20 30 40 50
                               60 70 80 90 100
Calculating graph cross linking in 332 Groups...
 % finished: 10 20 30 40 50
                               60 70 80 90 100
New number of ps-groups: 735
xsAnnotate has now 735 groups, instead of 332
Spectrum number increased after correlation done
Generating peak matrix for peak annotation!
Found and use user-defined ruleset!
Calculating possible adducts in 735 Groups...
 % finished: 10 20 30 40 50
                               60 70 80
Adduct/fragment annotation done
```

Because the parameter adductTable was not set in the peakAnnotation call, a warning was shown informing that the default MAIT table for positive polarisation mode was selected. The xsAnnotated object that contains all the information related to peaks, spectra and their annotation is stored in the MAIT object. It can be retrieved by typing:

#### > rawData(MAIT)

#### \$xsaFA

An "xsAnnotate" object!
With 735 groups (pseudospectra)
With 12 samples and 1261 peaks
Polarity mode is set to: positive
Using automatic sample selection
Annotated isotopes: 104
Annotated adducts & fragments: 72

# 4.4 Statistical Analysis

Memory usage: 4.43 MB

Following the first peak annotation stage, we want to know which features are different between classes. Consequently, we run the function spectralSigFeatures().

```
> MAIT<- spectralSigFeatures(MAIT.object = MAIT,pvalue=0.05,
         p.adj="none",scale=FALSE)
Skipping peak aggregation step...
> summary(MAIT)
A MAIT object built of 12 samples and 1261 peaks. No peak aggregation technique has been
56 of these peaks are statistically significant
The object contains 6 samples of class KO
The object contains 6 samples of class WT
Parameters of the analysis:
                                      Value
                                      "/home/biocbuild/bbs-3.20-bioc/R/site-library/faahKC
dataDir
snThres
                                      "2.12332257516562"
Sigma
mzSlices
                                      "0.3"
retcorrMethod
                                      "loess"
                                      "density"
groupMethod
                                      "3"
bwGroup
                                      "0.25"
mzWidGroup
filterMethod
                                      "centWave"
                                      "0.03"
rtStep
                                      "0"
nSlaves
                                      "MAIT Demo"
project
                                      "10"
ppm
                                      "0.5"
minfrac
fwhm
                                      "30"
                                      "gaussian"
family1
                                      "symmetric"
family2
                                      "0.2"
span
                                      "5"
centWave peakwidth1
                                      "20"
centWave peakwidth2
corrWithSamp
                                      "0.7"
                                      "0.75"
corrBetSamp
perfwhm
                                      "0.6"
                                      "6"
sigma
                                      "0.05"
peakAnnotation pvalue
                                      "TRUE"
calcIso
calcCiS
                                      "TRUE"
calcCaS
                                      "TRUE"
                                      "hcs"
graphMethod
```

"TRUE"

"None"

 ${\tt annotateAdducts}$ 

peakAggregation method

```
peakAggregation PCAscale "FALSE"
peakAggregation PCAcenter "FALSE"
peakAggregation scale "FALSE"
peakAggregation RemoveOnePeakSpectra "FALSE"
Welch pvalue "0.05"
Welch p.adj "none"
```

It is worth mentioning that by setting the scale parameter to TRUE, the data will be scaled to have unit variance. A summary of the statistically significant features is created and saved in a table called significantFeatures.csv (see Table ??). It is placed inside the Tables subfolder located in the project folder. This table shows characteristics of the statistically significant features, such as their P-value, the peak annotation or the expression of the peaks across samples. This table can be retrieved at any time from the MAIT-class objects by typing the instruction:

```
> signTable <- sigPeaksTable(MAIT.object = MAIT, printCSVfile = FALSE)
head(signTable)</pre>
```

The number of significant features can be retrieved from the MAIT-class object as follows:

> MAIT

A MAIT object built of 12 samples and 1261 peaks. No peak aggregation technique has been 56 of these peaks are statistically significant
The object contains 6 samples of class KO
The object contains 6 samples of class WT

#### 4.5 Statistical Plots

Out of 2,402 features, 106 were found to be statistically significant. At this point, several MAIT functions can be used to extract and visualise the results of the analysis. Functions plotBoxplot, plotHeatmap, plotPCA and plotPLS automatically generate boxplots, heat maps and PCA/PLS score plot files in the project folder when they are applied to a MAIT object (see Table ??).

```
> plotBoxplot(MAIT)
> plotHeatmap(MAIT)

> MAIT<-plotPCA(MAIT,plot3d=FALSE)
> MAIT<-plotPLS(MAIT,plot3d=FALSE)
> PLSmodel <- model(MAIT, type = "PLS")
> PCAmodel <- model(MAIT, type = "PCA")</pre>
```

> PLSmodel

## Partial Least Squares

```
12 samples56 predictors
```

2 classes: 'KO', 'WT'

No pre-processing

Resampling: Bootstrapped (25 reps)

Summary of sample sizes: 12, 12, 12, 12, 12, 12, ...

Resampling results across tuning parameters:

ncomp	${ t Accuracy}$	Kappa
1	1	1
2	1	1
3	1	1

Accuracy was used to select the optimal model using the largest value. The final value used for the model was ncomp = 1.

## > pcaScores(MAIT)

	PC1	PC2	PC3	PC4	PC5	PC6
[1,]	-8.103931	5.51253909	0.4560729	2.6304313	-0.001166111	-0.14024723
[2,]	-9.018332	-0.22222743	-3.7010516	-2.2074370	-1.944051226	0.65548080
[3,]	-5.511187	0.01698315	1.2817414	-2.5893508	3.427569605	-1.98128060
[4,]	-2.804242	-1.27431307	2.6576958	-0.8357823	0.453653596	3.36170251
[5,]	-2.964179	-4.19515061	0.7141187	2.3295230	-0.351542125	-1.13274686
[6,]	-3.879094	-2.88873255	0.1555561	1.3969922	-0.732631538	-0.30855890
[7,]	5.064055	1.00278256	0.6686535	0.3382944	0.542810579	1.22282372
[8,]	4.665604	0.94689331	1.7533831	-1.3360892	-2.774729880	-1.59218199
[9,]	5.978577	0.83212328	-0.4127802	-0.9160021	-0.177627419	0.26137007
[10,]	5.042800	0.59354948	0.9711509	-0.5728575	-1.200077540	-0.47892310
[11,]	5.587120	-0.03508374	-1.5934350	1.3148452	0.960806146	-0.08086651
[12,]	5.942808	-0.28936346	-2.9511057	0.4474328	1.796985912	0.21342810
	PC	C7 P(	28 PC	C9 PC	C10 PC1	.1 PC12
[1,]	-0.1863426	39 -0.389980 <i>6</i>	80 -0.216821	L3 -0.069358	321 -0.1049290	9.540979e-16
[2,]	-0.4250730	0.4423045	55 0.257178	39 0.232863	357 0.2261922	21 4.969983e-16
[3,]	0.1647694	13 0.2208554	12 0.137077	78 0.071010	0.0138375	50 -3.205769e-15
[4,]	-0.2998939	90 0.0330131	13 -0.998457	75 -0.092428	344 -0.2159734	7 -1.425943e-15
[5,]	-2.2938979	99 -0.2556343	39 0.262538	30 0.194420	0.1448242	26 1.089406e-15
[6,]	3.0277504	19 -0.3445307	75 0.400962	23 -0.437060	)48 -0.2988496	36 1.269818e-15
[7,]	0.0138996	0.9210798	33 1.924515	52 -0.342977	62 1.0650159	98 -2.185752e-16
[8,]	-0.2762524	1 0.7800297	73 -0.730497	79 -1.191974	189 0.0451882	25 4.943962e-16
[9,]	-0.6140281	l7 -1.3149737	76 1.119911	12 0.042059	945 -1.5495426	35 2.567391e-16
[10,]	0.7031050	9 -1.1230646	6 -0.433859	90 1.789984	190 0.7236246	50 5.846018e-16

```
[11,] 0.25951393 2.42089320 -0.6828096 0.73329844 -0.65932865 3.469447e-16 [12,] -0.07355034 -1.38999168 -1.0397381 -0.92983818 0.60994068 -3.053113e-16
```

The plotPCA and plotPLS functions produce MAIT objects with the corresponding PCA and PLS models saved inside. The models, loadings and scores can be retrieved from the MAIT objects by using the functions model, loadings and scores:

All the output figures are saved in their corresponding subfolders contained in the project folder. The names of the folders for the boxplots, heat maps and score plots are Boxplots, Heatmaps, PCA\_Scoreplots and PLS\_Scoreplots respectively. Inside the R session, the project folder is recovered by typing:

```
> resultsPath(MAIT)
```

[1] "/tmp/RtmpOdMmFg/Rbuild899c617d656f2/MAIT/vignettes/Results\_MAIT\_Demo"

## 4.6 Biotransformations

Before identifying the metabolites, peak annotation can be improved using the function Biotransformations to make interpreting the results easier. The MAIT package uses a default biotransformations table, but another table can be defined by the user and introduced by using the bioTable function input variable. The biotransformations table that MAIT uses is saved inside the file MAITtables.RData, under the name biotransformationsTable.

```
> Biotransformations(MAIT.object = MAIT, peakPrecision = 0.005)
```

```
% Annotation in progress: 10 20 30 40 50 60 70 80 90 100 A MAIT object built of 56 of these peaks are statistically significant. The object contains 6 samples of class KO. The object contains 6 samples of class WT.
```

Building a user-defined biotransformations table from the MAIT default table or adding a new biotransformation is straightforward. For example, let's say we want to add a new adduct called "custom" biotrans" whose mass loss is 105.

```
> data(MAITtables)
> myBiotransformation <- c ("custom_biotrans", 105.0)
> myBiotable<-biotransformationsTable
> myBiotable[,1]<-as.character(myBiotable[,1])</pre>
> myBiotable<-rbind(myBiotable,myBiotransformation)
> myBiotable[,1]<-as.factor(myBiotable[,1])</pre>
> tail(myBiotable)
                           NAME MASSDIFF
45
       glucuronide conjugation 176.0321
46 hydroxylation + glucuronide 192.027
                GSH conjugation 305.0682
47
48
    2x glucuronide conjugation 352.0642
49
                          [C13]
                                   1.0034
50
               custom_biotrans
                                      105
```

To build an entire new biotransformations table, you only need to follow the format of the biotransformations Table, which means writing the name of the biotransformations as factors in the NAME field of the data frame and their corresponding mass losses in the MASSDIFF field.

## 4.7 Metabolite Identification

Once the biotransformations annotation step is finished, the significant features have been enriched with a more specific annotation. The annotation procedure performed by the Biotransformations() function never replaces the peak annotations already done by other functions. MAIT considers the peak annotations to be complementary; therefore, when new annotations are detected, they are added to the current peak annotation and the identification function may be launched to identify the metabolites corresponding to the statistically significant features in the data.

```
> MAIT <- identifyMetabolites(MAIT.object = MAIT, peakTolerance = 0.005)
```

WARNING: No input database table was given. Selecting default MAIT database... Metabolite identification initiated

% Metabolite identification in progress: 10 20 30 40 50 60 70 80 90 100 Metabolite identification finished

By default, the function identifyMetabolites() looks for the peaks of the significant features in the MAIT default metabolite database. The input parameter peakTolerance defines the tolerance between the peak and a database compound to be considered a possible match. It is set to 0.005 mass/charge units by default. To check the results easily, function identifyMetabolites creates a table containing the significant feature characteristics and the possible metabolite identifications. Such a table is recovered from the MAIT-class object using the instruction:

- > metTable<-metaboliteTable(MAIT)
- > head(metTable)

	Query Ma	ass Databa	ase Mass	(neutral mass	s) rt	Isotop	ре	Adduct	Name
1	328	. 20		Unknov	n 56.32				Unknown
2	544	. 20		Unknov	n 56.62		[M+Na]+	521.202	Unknown
3	481	. 10		Unknov	n 61.42				Unknown
4	570	. 15		Unknov	n 53.81				Unknown
5	396	. 20		Unknov	n 64.39				Unknown
6	449	. 10		Unknov	n 54.87				Unknown
	spectra	${\tt Biofluid}$	ENTRY	p.adj		p F	Fisher.Test	Mean C	lass KO
1	4	unknown	unknown	0.0020024932	0.002002	24932	NA	49	9713.15
2	12	unknown	unknown	0.6602191037	0.660219	91037	NA	150	0448.04
3	429	unknown	unknown	0.0228596524	0.022859	96524	NA	32	2493.60
4	461	unknown	unknown	0.0219567137	0.021956	67137	NA	7:	1808.43
5	469	unknown	unknown	0.0178739635	0.017873	39635	NA	35	5365.04

```
6
          unknown unknown 0.0006846827 0.0006846827
                                                                  NA
                                                                          55495.48
 Mean Class WT Median Class KO Median Class WT KO WT
                                                             ko15
                                                                       ko16
1
       5636.569
                        42096.59
                                                         42881.87 88615.33
                                         4838.527
                                                       0
2
                                                       2 87514.52 15031.62
     117663.114
                       104227.29
                                       119004.818
                                                    3
3
       7120.502
                        36248.51
                                            0.000
                                                    3
                                                         15343.26 28757.71
4
      14649.083
                        78043.42
                                        13135.833
                                                    4
                                                             0.00 68659.67
5
       6391.812
                        30177.49
                                         6912.297
                                                    3
                                                       0
                                                         68932.53 49774.82
6
       3209.048
                        50442.30
                                                         57917.52 50601.15
                                         3153.240
                                                    6
                                                       0
                                                   wt15
       ko18
                ko19
                           ko21
                                       ko22
                                                             wt16
                                                                        wt18
   41311.31 38368.46
                       40095.58
                                  47006.340
                                               2120.425
                                                            0.000 11778.497
1
2 120940.07 70558.05 162788.17 445855.802 196595.300
                                                         6122.398 38846.610
   51219.32 50357.21
                       43739.32
                                   5544.795
                                             24105.695 18617.318
                                                                       0.000
4
  48512.19 87427.16 120206.09 106045.448
                                                  0.000 12884.211 13387.454
   29862.17 13673.40
                       19454.51
                                  30492.806
                                               4434.805
                                                         8914.240
                                                                    5780.494
   38995.65 45976.37
                      89198.74
                                  50283.450
                                               1360.024
                                                         2721.630
                                                                       0.000
       wt19
                   wt21
                              wt22
1 10243.441
               2888.511
                          6788.542
2 96823.161 141186.475 226404.742
3
      0.000
                  0.000
                             0.000
4 21632.388
                  0.000
                         39990.445
5
  8155.215
               3022.015
                          8044.100
  7280.855
              4306.930
                          3584.850
```

This table provides useful results about the analysis of the samples, such as the P-value of the statistical test, its adduct or isotope annotation and the name of any possible hit in the database. Note that if no metabolite has been found in the database for a certain feature, it is labelled as "unknown" in the table.

## 4.8 Validation

Finally, we will use the function Validation() to check the predictive value of the significant features. All the information related to the output of the Validation() function is saved in the project directory in a folder called "Validation". Two boxplots showing the overall and per class classification ratios are created, along with every confusion matrix corresponding to each iteration (see Table ??).

```
> MAIT <- Validation(Iterations = 20, trainSamples= 3,
+ MAIT.object = MAIT)

Iteration 1 done
Iteration 2 done
Iteration 3 done
Iteration 4 done
Iteration 5 done
Iteration 6 done
Iteration 7 done</pre>
```

Iteration 8 done
Iteration 9 done
Iteration 10 done
Iteration 11 done
Iteration 12 done
Iteration 13 done
Iteration 14 done
Iteration 15 done
Iteration 16 done
Iteration 17 done
Iteration 18 done
Iteration 19 done
Iteration 20 done

A summary of a MAIT object, which includes the overall classification values, can be accessed:

#### > summary(MAIT)

A MAIT object built of 12 samples and 1261 peaks. No peak aggregation technique has been 56 of these peaks are statistically significant. The object contains 6 samples of class KO

The object contains 6 samples of class WT

The Classification using 3 training samples and 20 Iterations gave the results:

 KNN PLSDA SVM

 mean
 1
 1
 1

 standard error
 0
 0
 0

#### Parameters of the analysis:

Value dataDir "/home/biocbuild/bbs-3.20-bioc/R/site-library/faahKC snThres Sigma "2.12332257516562" mzSlices "0.3" "loess" retcorrMethod groupMethod "density" "3" bwGroup "0.25" mzWidGroup "centWave"

filterMethod "centWave rtStep" "0.03" nSlaves "0"

project "MAIT\_Demo"

ppm "10" minfrac "0.5"

fwhm	"30"		
family1	"gaussian"		
family2	"symmetric"		
span	"0.2"		
centWave peakwidth1	"5"		
centWave peakwidth2	"20"		
corrWithSamp	"0.7"		
corrBetSamp	"0.75"		
perfwhm	"0.6"		
sigma	"6"		
peakAnnotation pvalue	"0.05"		
calcIso	"TRUE"		
calcCiS	"TRUE"		
calcCaS	"TRUE"		
graphMethod	"hcs"		
annotateAdducts	"TRUE"		
peakAggregation method	"None"		
peakAggregation PCAscale	"FALSE"		
peakAggregation PCAcenter	"FALSE"		
peakAggregation scale	"FALSE"		
<pre>peakAggregation RemoveOnePeakSpectra</pre>	"FALSE"		
Welch pvalue	"0.05"		
Welch p.adj	"none"		
peakTolerance	"0.005"		
polarity	"positive"		
Validation Iterations	"20"		
Validation trainSamples	"3"		
Validation PCAscale	"0"		
Validation PCAcenter	"1"		
Validation RemoveOnePeakSpectra	"0"		
Validation tuneSVM	"0"		
Validation scale	"1"		
PCA data logarithm	"FALSE"		
PCA data centered	"TRUE"		
PCA data scaled	"TRUE"		

It is also possible to gather the classification ratios per class, classifier used and iteration number by using the function classifRatioClasses():

## > classifRatioClasses(MAIT)

	KNN_Class_KO	PLSDA_Class_KO	SVM_Class_KO	KNN_Class_WT	PLSDA_Class_WT
[1,]	1	1	1	1	1
[2,]	1	1	1	1	1
[3,]	1	1	1	1	1
[4,]	1	1	1	1	1

```
[5,]
                      1
                                           1
                                                             1
                                                                               1
                                                                                                   1
 [6,]
                      1
                                           1
                                                             1
                                                                               1
                                                                                                   1
 [7,]
                      1
                                           1
                                                             1
                                                                               1
                                                                                                   1
 [8,]
                      1
                                                             1
                                           1
                                                                               1
                                                                                                   1
 [9,]
                      1
                                           1
                                                             1
                                                                                                   1
                                                                               1
[10,]
                      1
                                           1
                                                             1
                                                                                                   1
                                                                               1
[11,]
                      1
                                           1
                                                             1
                                                                               1
                                                                                                   1
[12,]
                      1
                                           1
                                                             1
                                                                                                   1
                                                                               1
[13,]
                      1
                                           1
                                                             1
                                                                               1
                                                                                                   1
[14,]
                      1
                                           1
                                                             1
                                                                               1
                                                                                                   1
[15,]
                      1
                                           1
                                                             1
                                                                               1
                                                                                                   1
[16,]
                      1
                                           1
                                                             1
                                                                               1
                                                                                                   1
[17,]
                      1
                                           1
                                                             1
                                                                                                   1
                                                                               1
[18,]
                      1
                                           1
                                                             1
                                                                               1
                                                                                                   1
[19,]
                      1
                                           1
                                                             1
                                                                               1
                                                                                                   1
[20,]
                      1
                                           1
                                                             1
                                                                               1
                                                                                                   1
       SVM_Class_WT
 [1,]
                      1
 [2,]
                      1
                      1
 [3,]
                      1
 [4,]
 [5,]
                      1
 [6,]
                      1
 [7,]
                      1
 [8,]
                      1
 [9,]
                      1
                      1
[10,]
[11,]
                      1
[12,]
                      1
[13,]
                      1
[14,]
                      1
[15,]
                      1
[16,]
                      1
                      1
[17,]
                      1
[18,]
                      1
[19,]
[20,]
                      1
```

The classification ratios are 100%; the set of significant features separates the samples belonging to these classes.

# 4.9 Using External Peak Data

Taking advantage of the modularised design of MAIT, it is possible to use the function MAIT-builder to import peak data and analyse it using the MAIT statistical functions. As stated in

section ??, there are certain arguments that should be provided depending on which function is wanted to be launched. In this section we will show an example of this data importation procedure using the same data that we have been using in the tutorial so far. Let's say we have a peak table recorded in positive polarisation mode with the peak masses and retention time values such as:

```
> peaks <- scores(MAIT)
> masses <- getPeaklist(MAIT)$mz
> rt <- getPeaklist(MAIT)$rt/60</pre>
```

We want to perform an annotation stage and metabolite identification on these data. To that end, we can launch the function MAITbuilder to build a MAIT-class object with the data in the table:

We have selected the option spectraEstimation as TRUE because we do not know the grouping of the peaks into spectra. As we want to annotate and identify all the peaks in the data frame, we set the flag significantFeatures to TRUE. At this point, we can launch the Biotransformations function:

```
> importMAIT <- Biotransformations(MAIT.object = importMAIT,
+ adductAnnotation = TRUE,
+ peakPrecision = 0.005, adductTable = NULL)</pre>
```

Set adductTable equal to negAdducts to use the default MAIT table for negative polarity

```
% Annotation in progress: 0 10 20 30 40 50 60 70 80 90 100
```

We set the adductAnnotation flag to TRUE as we want to perform an adduct annotation step. The parameter adductTable set to NULL implies that a positive polarisation adduct annotation stage will be performed. To run a negative annotation, the argument should be set to negAdducts. The metabolite identification stage is launched as in the previous case:

WARNING: No input database table was given. Selecting default MAIT database... Metabolite identification initiated

% Metabolite identification in progress: 0 10 20 30 40 50 60 70 80 90 100 Metabolite identification finished